

Selection Manual

EFF Dental Screws



MULTI PLATFORM SYSTEMS



Let's understand a
little more about
screws?



Hex Screw



Dynamic
Screw



EGM
screw



Grip Screw



Torx screw



Square
Screw



And these torque screw
code tips!
Come see!



Legend:



Systems ready for ~~del~~ **del** systems

EFF 2101...
10Ncm

EFF 2102...
20Ncm

EFF 2103...
20Ncm


EFF 2104...
32Ncm

EFF 2105...
32Ncm




Solutions for all systems


With different connections, EFF Dental screws meet all systems on the national and international market.

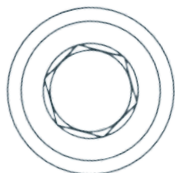


Hexagon

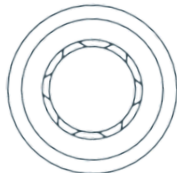


GRIP

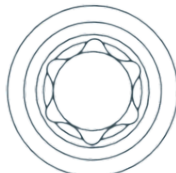


Trox

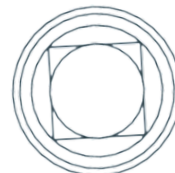

Square

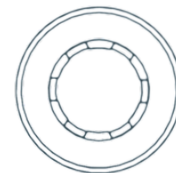

EGM

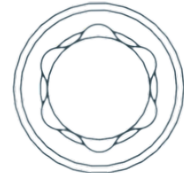

360°



EFF[®]
DENTAL COMPONENTES





Nobel Biocare® 3.0	K3.1	pg.5
Nobel Biocare® 3.5 NP	K3.2	pg.6
Nobel Biocare® 4.3 RP	K3.3	pg. 7
Nobel Biocare® 5.5 WP	K3.5	pg. 8
Straumann® 4.8 RN	J3.1	pg. 9
Straumann® 6.5 WN	J3.2	pg.10
Straumann® 2.9 SC	3.5	pg. 11
Straumann® 3.3 NC	J3.3	pg. 12
Straumann® 4.1 RC	J3.4	pg. 13
Straumann® 3.5/3.75/4.0/4.5 RB	J3.6	pg.14
Straumann® 3.5/3.75/4.0/4.55.0/5.5/6.5 WB	J3.7	pg.15
Straumann® 3.5 NT	J3.9	pg. 16
Straumann® 4.8 RT	J3.10	pg. 16
Straumann® 6.5 WT	J3.11	pg. 16
Pilar Syno RN	SYRN	pg. 17
Pilar Syno WN	SYWN	pg. 18
Pilar MBase 3.5 4.5	MB	pg. 19
Pilar SRA 3.5 4.6	SRA	pg. 20
Neodent® 3.5 3.7 4.0 4.3 5.0 CM	A3.1	pg. 21
Neodent® 3.0 EGM	A3.2	pg. 22
Pilar EGM	PGM	pg. 23
Micro Pilar EGM	MPGM	pg. 24
Pilar GT	MPGT	pg. 25
SIN® Strong® 3.5/3.8/4.5/5.5	C3.1	pg. 26
MIS® CC NP C1 Narrow	T3.1	pg. 27
MIS® CC NP V3 Narrow	T3.2	pg. 28
MIS® CC SP C1/V3 Standard	T3.3	pg. 29
Plenum® 3.0 Slim	Y3.1	pg. 30
Plenum® 3.5 REGULAR	Y3.2	pg. 31
Bionnovation 3.5 Biomorse	U3.1	pg. 32
Implacil® Due Cone	D3.1	pg. 33
Signo Vincés® Duocon	Q3.1	pg. 34
Astra Tech® 4.5 5.0 Ossoespeed	Q3.2	pg. 35
Astra Tech® 3.6 EV	Q3.3	pg. 36
Astra Tech® 4.2 EV	Q3.4	pg. 37
Astra Tech® 4.8 EV	Q3.5	pg. 38
Astra Tech® 5.4 EV	Q3.6	pg. 39



Nobel Biocare® 3.5 NP	K2.1	pg. 41
Nobel Biocare® 4.3 RP	K2.2	pg. 42
Nobel Biocare® 5.0 WP	K2.3	pg. 43
Implacil® 3.5	D2.1	pg. 44
Implacil® 4.0	D2.2	pg. 45
Implacil® 5.0	D2.3	pg. 46
Frialit® 3.8 SIN® 3.8	O2.2 B2.1	pg. 47
Frialit® 4.5 SIN® 4.5	O2.3 B2.2	pg. 48
Frialit® 5.5 SIN® 5.5	O2.4 B2.3	pg. 49
Dentoflex® 3.25	G2.1	pg. 50
Dentoflex® 3.5	G2.2	pg. 51
Dentoflex® 3.75	G2.3	pg. 52
Camlog® 3.3	P2.1	pg. 53
Camlog® 3.8	P2.2	pg. 54
Camlog® 4.3	P2.3	pg.55
Camlog® 5.0	P2.4	pg. 56
Camlog® 6.0	P2.5	pg. 57
Baumer® 4.0	X2.1	pg. 58
Conexão® 4.1/4.3/5.0/6.0	C2.1 C2.2 C2.3	pg. 59
Titaniumfix® 3.75/3.0	E2.1	pg. 60
Biomet 3i™ 3.4	L2.1	pg. 61
Biomet 3i™ 4.1	L2.2	pg. 62
Biomet 3i™ 5.0	L2.3	pg. 63
Neodent® 4.3 5.0	A2.1	pg. 64



Branemark System® 3.5	K1.1	pg. 66
Branemark System® 4.1	K1.2	pg. 67
Branemark System® 3.5	K1.3	pg. 68
Neodent® 3.3 Conexão® 3.3	A1.1 C1.1	pg. 69
Neodent® 4.3	A1.3	pg. 70
Neodent® 5.0 Conexão® 5.0	A1.4 C1.3	pg. 71
Biomet 3i™ 3.4 SIN® 3.4	L1.1?B1.1	pg. 72
Biomet 3i™ 5.0 SIN® 5.0	L1.3 B1.3	pg. 73
implacil® 3.3	D1.1	pg. 74
Titaniumfix® 3.3	E1.1	pg. 75



INTERMEDIATE PILLAR

Multi-unit - Regular Platform	M.P.	pg. 76
Multi-unit - Long Platform	M.P.L1	pg.77
Multi-unit - Long Platform - NEO	M.P.L2	pg. 78
Esthetione	EONE	pg. 79
Estheticone	EST	pg. 80



Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2101.06A



Grip
EFF2101.06C



10Ncm



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2101.06A



Grip
EFF2101.06C



10Ncm



Prov. EVOLUTION



Hexagonal
EFF2101.06A



Grip
EFF2101.06C



10Ncm



Link to CAD/CAM TiBase



Hexagonal
EFF2101.06A



Grip
EFF2101.06C



10Ncm



Dynamic Link



Hexagonal
EFF2101.06A



Grip
EFF2101.06C



10Ncm



Interface



Grip
EFF2101.07C



10Ncm



K3.1

NOBEL BIOCARE® 3.0

CONE MORSE





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal!
Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2102.02A



Grip
EFF2102.02C



Dynamic
EFF2102.02F



Ucla Plastic
Surgery with CrCo



Hexagonal
EFF2102.02A



Grip
EFF2102.02C



Dynamic
EFF2102.02F



K3.2

**NOBEL BIOCARE®
3.5|NP**

CONE MORSE



**Prov.
EVOLUTION**



Hexagonal
EFF2102.02A



Grip EFF2102.02C



Dynamic
EFF2102.02F



**Link to
CAD/CAM
TiBase**



Hexagonal
EFF2102.02A



Grip
EFF2102.02C



Dynamic
EFF2102.02F



**Dynamic
Link**



Hexagonal
EFF2102.02A



Grip
EFF2102.02C



Dynamic
EFF2102.02F



Interface



Grip
EFF2102.08C





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2104.15A



Grip
EFF2104.15C



Dynamic
EFF2104.15F

32Ncm



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2104.15A



Grip
EFF2104.15C



Dynamic
EFF2104.15F

32Ncm



K3.3

NOBEL BIOCARE® 4.3|RP

CONE MORSE



Prov. EVOLUTION



Hexagonal
EFF2104.15A



Grip
EFF2104.15C



Dynamic
EFF2104.15F

32Ncm



Link to CAD/CAM TiBase



Hexagonal
EFF2104.15A



Grip
EFF2104.15C



Dynamic
EFF2104.15F

32Ncm



Dynamic Link



Hexagonal
EFF2104.15A



Grip
EFF2104.15C



Dynamic
EFF2104.15F

32Ncm




Interface



Grip
EFF2104.13C

32Ncm

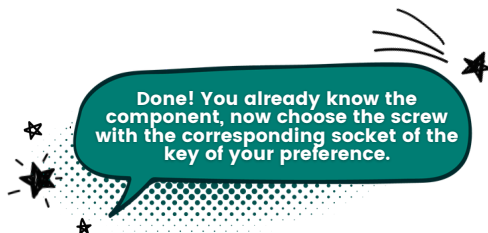

K3.5

NOBEL BIOCARE®
5.5|WP

CONE MORSE



Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.15A



Grip
EFF2104.15C



Dynamic
EFF2104.15F



Dynamic
Link




Hexagonal
EFF2104.15A



Grip
EFF2104.15C




Dynamic
EFF2104.15F


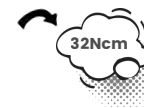


Interface



Grip
EFF2104.13C

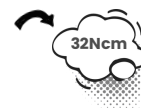


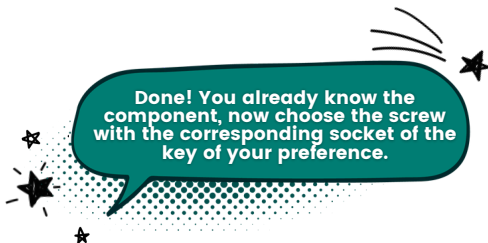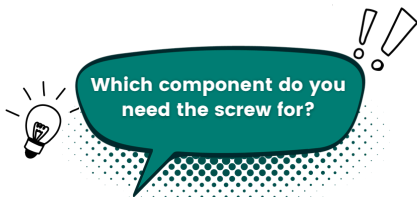

Temporary Abutment



Torx EFF2104.11D
S/AR EFF2104.17D
C/AR

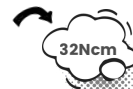

Ucla Plastic Surgery with CrCo



Torx EFF2104.11D
S/AR EFF2104.17D
C/AR



Prov. EVOLUTION



Torx EFF2104.17D
C/AR



Link to CAD/CAM TiBase



Torx EFF2104.17D
C/AR



Link to CAD/CAM Perfect Base





Torx EFF2104.19D



J3.1

STRAUMANN®
4.8|RN
CONE MORSE





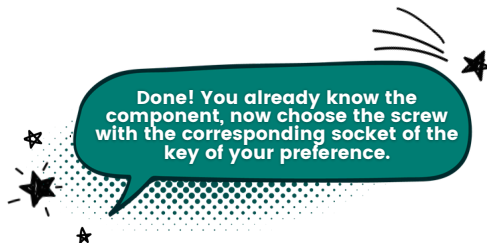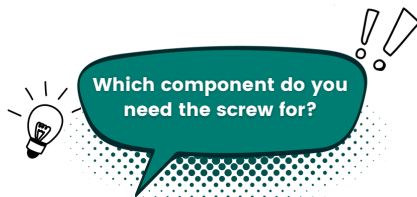

J3.2

STRAUMANN®
6.5|WN
CONE MORSE



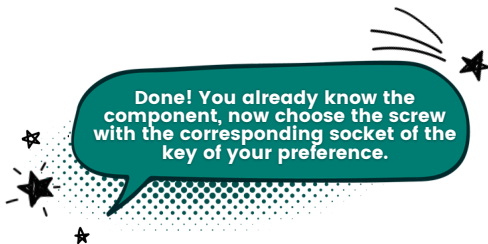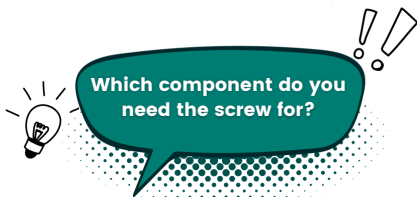

System Multi Plattform

J3.5

STRAUMANN®
29|SC
CONE MORSE



Temporary
Abutment



Ucla Plastic
Surgery with CrCo



Prov.
EVOLUTION



Link to
CAD/CAM
TiBase



Interface



Torx EFF2101.13D



Torx EFF2101.13D



Torx EFF2101.13D



Torx EFF2101.13D



Torx EFF2101.15D





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Torx EFF2102.11D



Dynamic EFF2102.11F



20Ncm



Ucla Plastic Surgery with CrCo



Torx EFF2102.11D



Dynamic EFF2102.11F



20Ncm



Prov. EVOLUTION



Torx EFF2102.11D



Dynamic EFF2102.11F



20Ncm



Link to CAD/CAM TiBase



Torx EFF2102.11D



Dynamic EFF2102.11F



20Ncm



Link to CAD/CAM Perfect Base



Torx EFF2102.11D



Dynamic EFF2102.11F



20Ncm



Interface



Torx EFF2102.05D



20Ncm



J3.3

STRAUMANN®
3.3|NC

CONE MORSE





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal!
Here you check the correct torque of the screw



Temporary Abutment



Torx EFF2102.20D



Dynamic EFF2102.11F



Ucla Plastic Surgery with CrCo



Torx EFF2102.20D



Dynamic EFF2102.11F



Prov. EVOLUTION



Torx EFF2102.20D



Dynamic EFF2102.11F



Link to CAD/CAM TiBase



Torx EFF2102.20D



Dynamic EFF2102.11F



Link to CAD/CAM Perfect Base



Torx EFF2102.20D



Dynamic EFF2102.11F



Interface



Torx EFF2102.05D





J3.4

STRAUMANN®
4i|RC

CONE MORSE



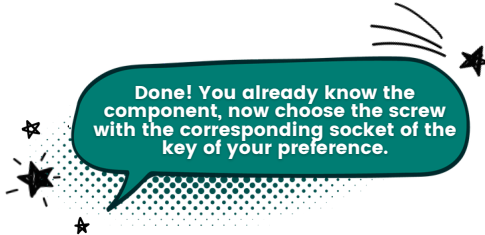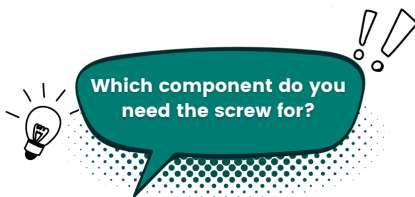

Temporary Abutment



Torx EFF2102.45D



Ucla Plastic Surgery with CrCo



Torx EFF2102.45D



Prov. EVOLUTION



Torx EFF2102.45D



Link to CAD/CAM TiBase



Torx EFF2102.45D



Link to CAD/CAM Perfect Base



Torx EFF2102.45D



Interface





Torx EFF2102.47D



J3.6

STRAUMANN®
3.5/3.75/4.0/4.5|RB
CONE MORSE



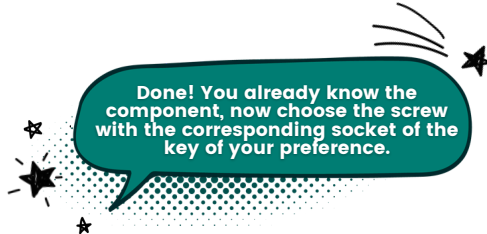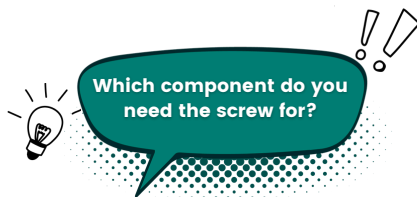

Temporary Abutment



Torx EFF2102.45D



Ucla Plastic Surgery with CrCo



Torx EFF2102.45D



Prov. EVOLUTION



Torx EFF2102.45D



Link to CAD/CAM TiBase



Torx EFF2102.45D



Link to CAD/CAM Perfect Base



Torx EFF2102.45D



Interface




Torx EFF2102.47D



J3.7

STRAUMANN®
3.5/3.75/4.0/4.5
5.0/5.5/6.5|WB
CONE MORSE



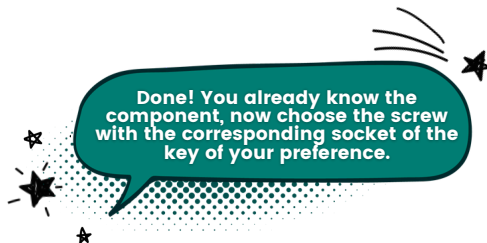


SMP
Sistema Multi-Platform
J3.9
STRAUMANN®
3.5|NT
CONE MORSE



Link to CAD/CAM
TiBase



Torx EFF2102.82D




SMP
Sistema Multi-Platform
J3.10
STRAUMANN®
4.8|RT
CONE MORSE



Link to CAD/CAM
TiBase



Torx EFF2102.82D




SMP
Sistema Multi-Platform
J3.11
STRAUMANN®
6.5|WT
CONE MORSE





Link to CAD/CAM
TiBase



Torx EFF2102.82D



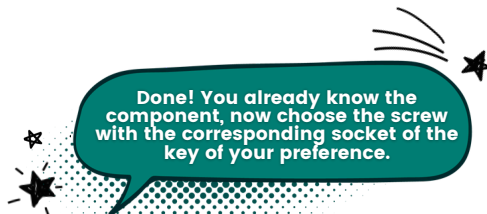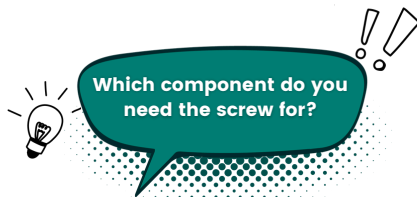

Plastic Coping



Torx EFF2104.10D



SYRN

PILAR SYNO RN

CONE MORSE



CrCo Base Plastic Coping



Torx EFF2104.10D





Link to CAD/CAM TiBase



Torx EFF2104.10D



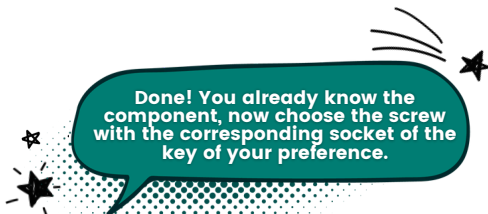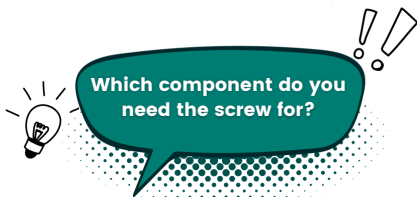

Plastic Coping



Torx EFF2104.10D



32Ncm



SYWN
PILAR SYNO WN

CONE MORSE



CrCo Base
Plastic Coping



Torx EFF2104.10D



32Ncm






Link to
CAD/CAM
TiBase



Torx EFF2104.10D



32Ncm

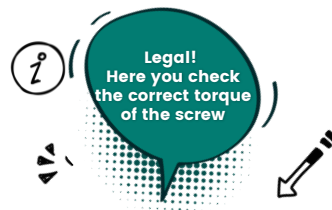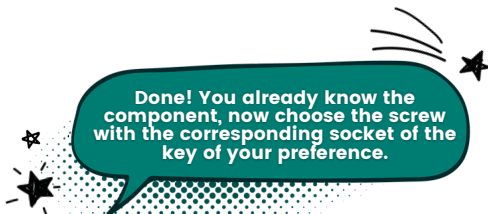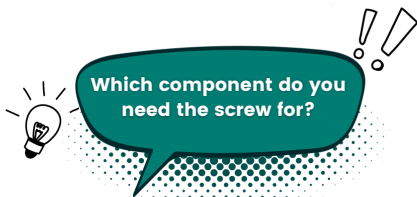

Plastic Coping



Torx EFF2101.08D



MB
MBASE PILLAR
3.5/4.5
CONE MORSE



CrCo Base Plastic Coping



Torx EFF2101.08D





Link to
CAD/CAM
TiBase



Torx EFF2101.08D




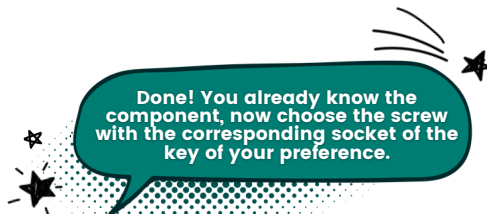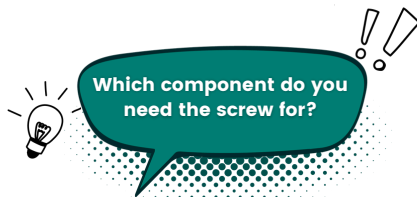

Plastic Coping



Torx EFF2101.11D

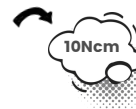

10Ncm



CrCo Base Plastic Coping



Torx EFF2101.11D



10Ncm



Link to CAD/CAM TiBase



Torx EFF2101.11D



10Ncm





SRA

PILAR SRA
3.5/4.6

CONE MORSE



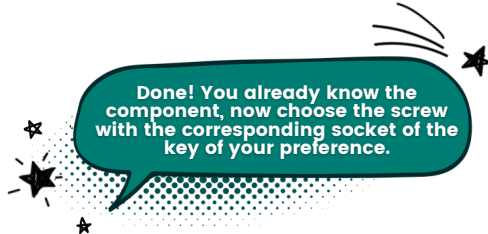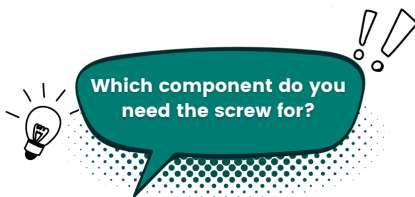

Temporary Abutment



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)
EFF2103.22A.3
(for h3.50/4.50mm)



Dynamic
EFF2103.22F.2
(for h0.80/1.50/2.50mm)
EFF2103.22F.3
(for h3.50/4.50mm)



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)
EFF2103.22A.3
(for h3.50/4.50mm)



Dynamic
EFF2103.22F.2
(for h0.80/1.50/2.50mm)
EFF2103.22F.3
(for h3.50/4.50mm)



Prov. EVOLUTION



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)
EFF2103.22A.3
(for h3.50/4.50mm)



Dynamic
EFF2103.22F.2
(for h0.80/1.50/2.50mm)
EFF2103.22F.3
(for h3.50/4.50mm)



Link to CAD/CAM TiBase



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)
EFF2103.22A.3
(for h3.50/4.50mm)



Dynamic
EFF2103.22F.2
(for h0.80/1.50/2.50mm)
EFF2103.22F.3
(for h3.50/4.50mm)



Dynamic Link



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)
EFF2103.22A.3
(for h3.50/4.50mm)



Dynamic
EFF2103.22F.2
(for h0.80/1.50/2.50mm)
EFF2103.22F.3
(for h3.50/4.50mm)



Interface



Hexagonal
EFF2103.09A



Sistema Multi Piattaforma

A3.1

NEODENT®
3.5/3.7/4.0/4.3/5.0
CM
CONE MORSE



Which component do you need the screw for?

Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.

Legal! Here you check the correct torque of the screw



Temporary Abutment



EGM EFF2102.32E



Ucla Plastic Surgery with CrCo



EGM EFF 2102 . 34E



Dynamic EFF 2102 . 34F



Prov. EVOLUTION



EGM EFF2102.32E



Cementable trunnion



EGM EFF2102.32E



Link to CAD/CAM TiBase



EGM EFF2102.32E



Dynamic Link



EGM EFF2102.32E



Dynamic EFF2102.34F



Interface



EGM EFF2102.41E



20Ncm



20Ncm



20Ncm



20Ncm



20Ncm



20Ncm





20Ncm



A3.2

NEODENT®
3.0 EGM
CONE MORSE




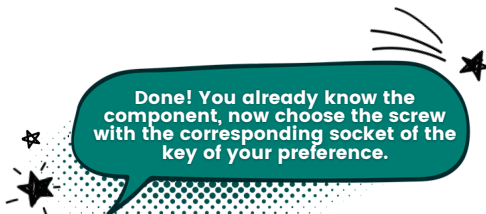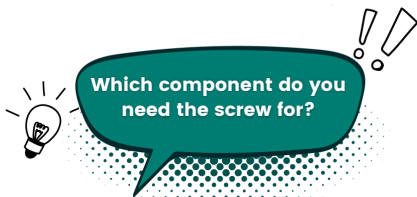

Temporary Abutment



EGM EFF2104.24E

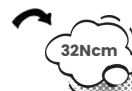

CrCo Base Plastic Coping




EGM EFF2104.24E



Link to CAD/CAM TiBase





EGM EFF2104.24E

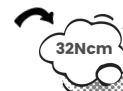

PGM
PILAR EGM

CONE MORSE



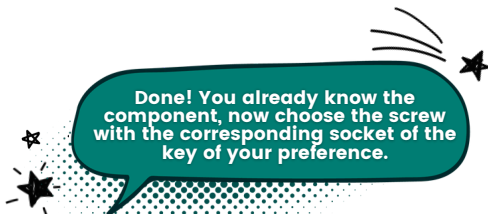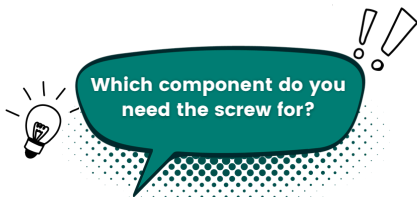

Temporary Abutment



Hexagonal
EFF2101.12A



CrCo Base
Plastic Coping



Hexagonal
EFF2101.12A



Link to
CAD/CAM
TiBase



Hexagonal
EFF2101.12A





MPGM

MICRO PILAR EGM

MICRO ABUTMENT



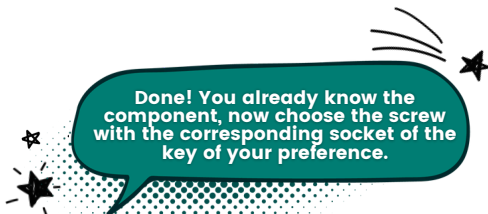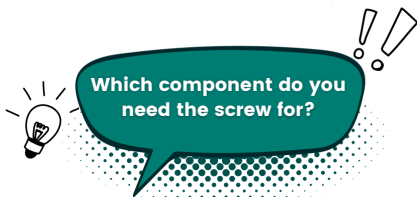

Temporary Abutment



Hexagonal
EFF2104.22A



CrCo Base
Plastic Coping



Hexagonal
EFF2104.22A



Link to
CAD/CAM
TiBase





Hexagonal
EFF2104.22A



MPGT
ABUTMENT GT

CONE MORSE



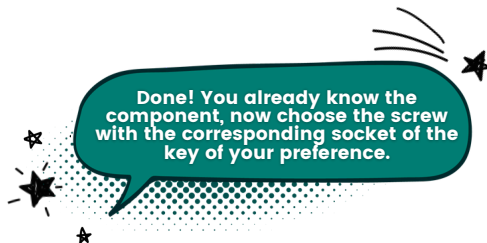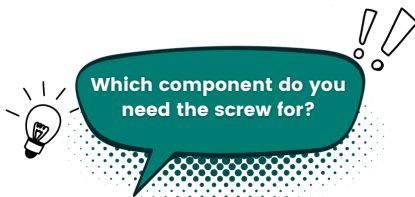

C3.1

SIN® STRONG®
3.5/3.8/4.5/5.5
CONE MORSE



Temporary
Abutment



Ucla Plastic
Surgery with CrCo



Temporary
Abutment



Link to
CAD/CAM
TiBase



Interface



Hexagonal
EFF2102.24A



Hexagonal
EFF2102.24A



Hexagonal
EFF2102.24A





Hexagonal
EFF2102.24A



Hexagonal
EFF2102.18A



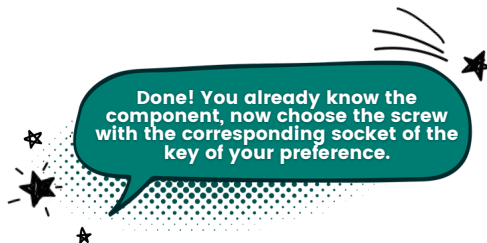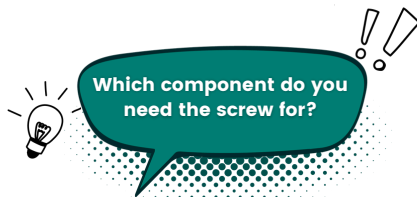

T3.1

MIS® CC NP C1
NARROW
CONE MORSE



Temporary
Abutment



Hexagonal 1.3mm
EFF2102.55A



Dynamic
EFF2102.55F



Ucla Plastic
Surgery with CrCo



Hexagonal 1.3mm
EFF2102.55A



Dynamic
EFF2102.55F



Prov.
EVOLUTION



Hexagonal 1.3mm
EFF2102.55A



Dynamic
EFF2102.55F



Link to
CAD/CAM
TiBase



Hexagonal 1.3mm
EFF2102.55A



Dynamic
EFF2102.55F





Interface



Hexagonal 1.3mm
EFF2102.54A



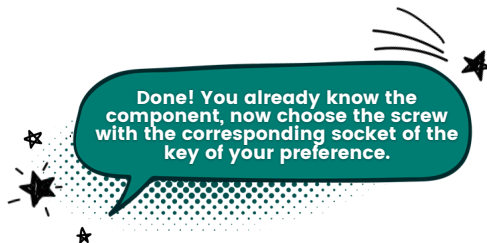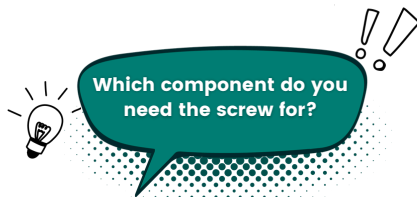

T3.2

MIS® CC NP V3
NARROW
CONE MORSE



Temporary
Abutment



Ucla Plastic
Surgery with CrCo



Prov.
EVOLUTION



Link to
CAD/CAM
TiBase



Interface



Hexagonal 1.3mm
EFF2102.55A



Dynamic
EFF2102.55F



Hexagonal 1.3mm
EFF2102.55A



Dynamic
EFF2102.55F



Hexagonal 1.3mm
EFF2102.55A



Dynamic
EFF2102.55F



Hexagonal 1.3mm
EFF2102.55A



Dynamic
EFF2102.55F



Hexagonal 1.3mm
EFF2102.54A





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Hexagonal 1.3mm
EFF2103.14A



Dynamic
EFF2103.14F



Ucla Plastic
Surgery with CrCo



Hexagonal 1.3mm
EFF2103.14A



Dynamic
EFF2103.14F



Prov.
EVOLUTION



Hexagonal 1.3mm
EFF2103.14A



Dynamic
EFF2103.14F



Link to
CAD/CAM
TiBase



Hexagonal 1.3mm
EFF2103.14A



Dynamic
EFF2103.14F



Dynamic
Link



Hexagonal 1.3mm
EFF2103.14A



Dynamic
EFF2103.14F



Interface



Hexagonal 1.3mm
EFF2103.42A





T3.3

MIS® CC SP C1/V3
STANDARD

CONE MORSE



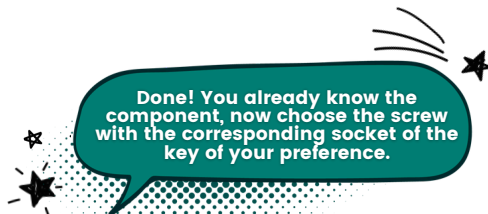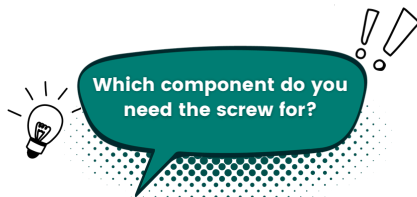

Y3.1

PLENUM®
3.0|SLIM

CONE MORSE



Prov.
EVOLUTION



Hexagonal
EFF2101.65A.1





Link to
CAD/CAM
TiBase



Hexagonal
EFF2101.65A.1
(1.0|2.0|3.0mm tapes)
EFF2101.65A.2
(tapes 4.0|5.0|6.0mm)



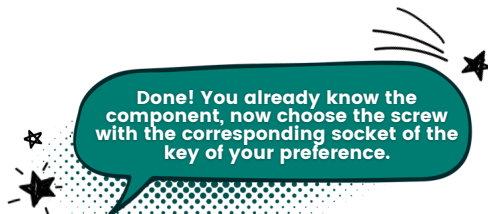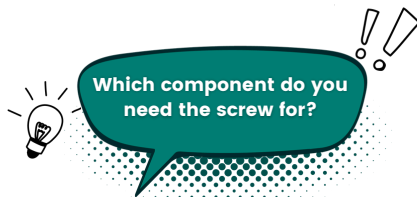

Y3.2

PLENUM®
3.5|REGULAR

CONE MORSE



Prov.
EVOLUTION



Hexagonal
EFF2102.66A.1





Link to
CAD/CAM
TiBase



Hexagonal
EFF2102.66A.1
(1.0|2.0|3.0mm tapes)
EFF2102.66A.2
(tapes 4.0|5.0|6.0mm)



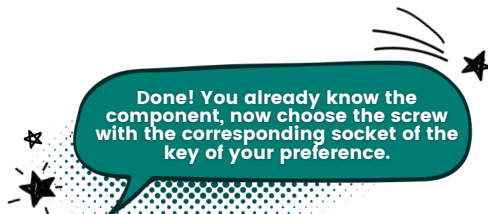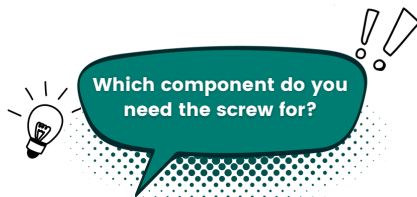

Prov.
EVOLUTION



Hexagonal
EFF2102.78A



U3.1

BIONNOVATION 3.5|BIOMORSE

CONE MORSE



Link to
CAD/CAM
TiBase



Hexagonal
EFF2102.78A





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal!
Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)
EFF2103.22A.3
(for h3.50/4.50mm)



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)
EFF2103.22A.3
(for h3.50/4.50mm)



Prov. EVOLUTION



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)
EFF2103.22A.3
(for h3.50/4.50mm)



Link to CAD/CAM TiBase



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)



Dynamic Link



Hexagonal
EFF2103.22A.2
(for h0.80/1.50/2.50mm)



Interface





Hexagonal
EFF2103.09A



D3.1
IMPLACIL®
TWO CONES
CONE MORSE



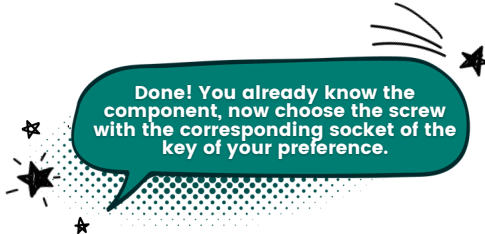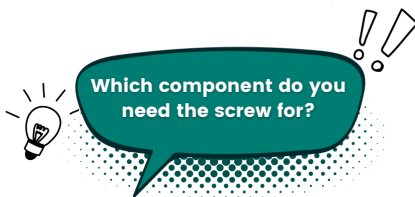

Temporary Abutment



**Hexagonal
EFF2102.06A**



20Ncm



Ucla Plastic Surgery with CrCo



**Hexagonal
EFF2102.06A**



20Ncm



**Prov.
EVOLUTION**



**Hexagonal
EFF2102.06A**



20Ncm



**Link to
CAD/CAM
TiBase**



**Hexagonal
EFF2102.06A**



20Ncm



**Dynamic
Link**



**Hexagonal
EFF2102.06A**



20Ncm



Interface



**Hexagonal
EFF2102.15A**





20Ncm



Q3.1

**VINCES® SIGN
DUOCON
CONE MORSE**



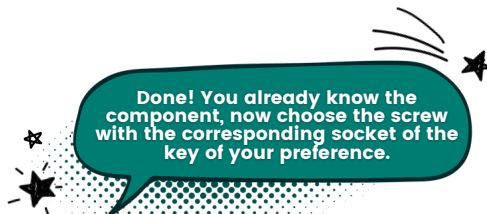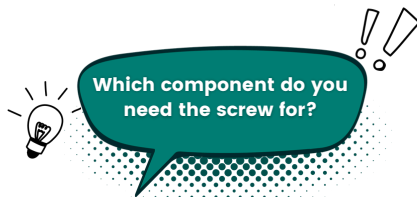

Q3.2

ASTRA TECH®
4.5|5.0 OSSOESPEED

CONE MORSE



Temporary
Abutment



Hexagonal
EFF2104.02A





Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.02A



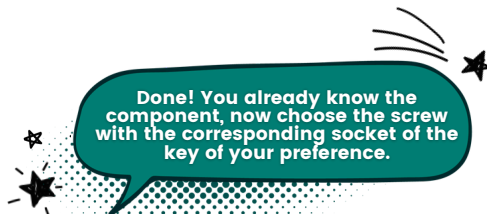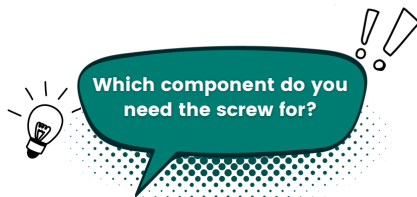

Q3.3

ASTRA TECH®
3.6|EV

CONE MORSE



Temporary
Abutment



Hexagonal
EFF2102.29A



20Ncm





Link to
CAD/CAM
TiBase



Hexagonal
EFF2102.29A



20Ncm

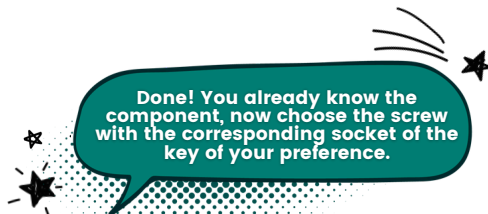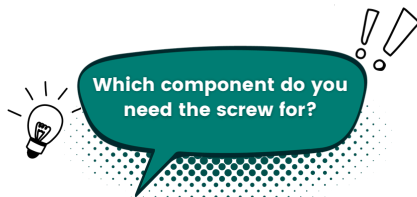

Q3.4

ASTRA TECH®
4.2|EV

CONE MORSE



Temporary
Abutment



Hexagonal
EFF2103.27A





Link to
CAD/CAM
TiBase



Hexagonal
EFF2103.27A



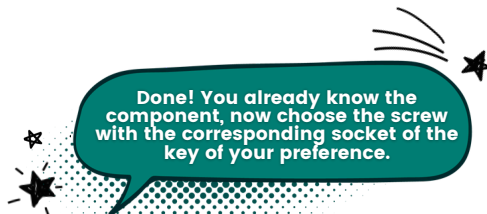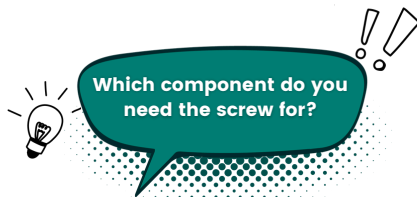

Q3.5

ASTRA TECH®
4.8|EV

CONE MORSE



Temporary
Abutment



Hexagonal
EFF2104.21A





Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.21A



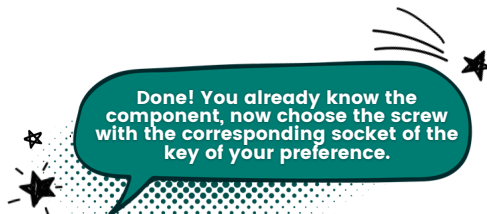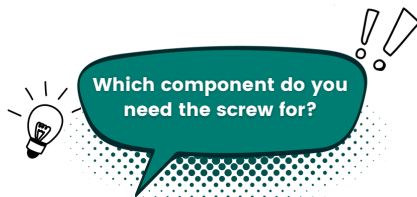

Q3.6

ASTRA TECH®
5.4|EV

CONE MORSE



Temporary
Abutment



Hexagonal
EFF2104.21A



Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.21A





3 REASONS



for you to start using
the **Driver and the Screw 360°**



1

Elimination of the
vestibular canal;


2

Tilt up to 12° at different
angles;



3

100% improvement in
dental aesthetics;

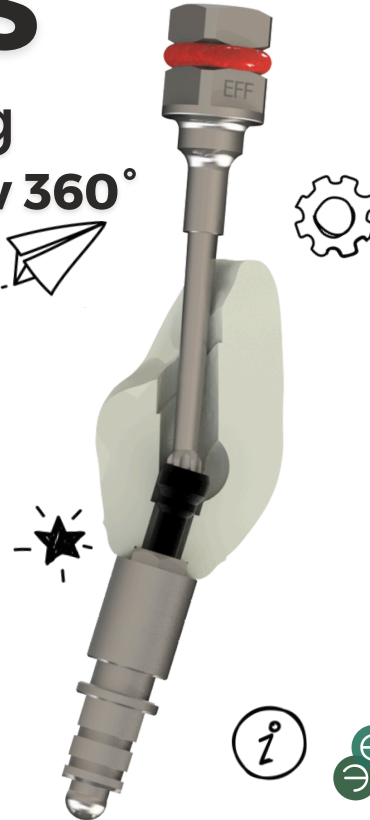

EEE[®]
DENTAL COMPONENTES





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2103.03A



Grip
EFF2103.03C



Dynamic
EFF2103.03F



Square
EFF2103.03B



Sistema Multi Plataforma

K2.1



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2103.03A



Grip
EFF2103.03C



Dynamic
EFF2103.03F



Square
EFF2103.03B



**NOBEL BIOCARE®
3.5|NP**

INTERNAL CONNECTION



PREMIUM



Prov. EVOLUTION



Hexagonal
EFF2103.03A



Grip
EFF2103.03C



Dynamic
EFF2103.03F



Square
EFF2103.03B



Interface



Grip
EFF2103.04C





Which component do you need the screw for?




Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal!
Here you check the correct torque of the screw



Temporary Abutment

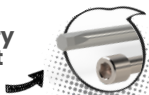

Hexagonal
EFF2104.02A




Grip
EFF2104.02C



Dynamic
EFF2104.02F



Square
EFF2104.02B

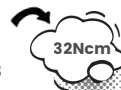

K2.2

**NOBEL BIOCARE®
43|RP**

INTERNAL CONNECTION



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2104.02A



Grip
EFF2104.02C



Dynamic
EFF2104.02F



Square
EFF2104.02B



Prov. EVOLUTION



Hexagonal
EFF2104.02A




Grip
EFF2104.02C



Dynamic
EFF2104.02F



Square
EFF2104.02B

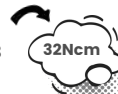

Interface



Grip
EFF2104.04C





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal!
Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2104.02A



Grip
EFF2104.02C



Dynamic
EFF2104.02F



Square
EFF2104.02B



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2104.02A



Grip
EFF2104.02C



Dynamic
EFF2104.02F



Square
EFF2104.02B



K2.3

NOBEL BIOCARE®
5.0|WP

INTERNAL CONNECTION



Prov. EVOLUTION



Hexagonal
EFF2104.02A



Grip
EFF2104.02C



Dynamic
EFF2104.02F



Square
EFF2104.02B



Interface



Grip
EFF2104.04C





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2103.02A



Square
EFF2103.02BB



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2103.02A



Square
EFF2103.02BB



Prov. EVOLUTION



Hexagonal
EFF2103.02A



Square
EFF2103.02BB



Link to CAD/CAM TiBase



Hexagonal
EFF2103.02A



Square
EFF2103.02BB



Dynamic Link



Hexagonal
EFF2103.02A



Square
EFF2103.02BB



Interface



Hexagonal
EFF2103.08A



D2.1

IMPLACIL® 3.5

INTERNAL CONNECTION





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw




Temporary Abutment



Hexagonal
EFF2104.01A



Square
EFF2104.01B

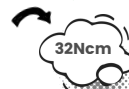


Ucla Plastic Surgery with CrCo



Hexagonal
EFF2104.01A



Square
EFF2104.01B

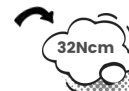


Prov. EVOLUTION



Hexagonal
EFF2104.01A



Square
EFF2104.01B

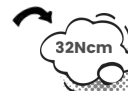


Link to CAD/CAM TiBase



Hexagonal
EFF2104.01A



Square
EFF2104.01B

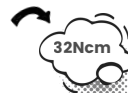

Dynamic Link



Hexagonal
EFF2104.01A



Square
EFF2104.01B



Interface



Hexagonal
EFF2104.03A



D2.2
IMPLACIL®
4.0

INTERNAL
CONNECTION





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2104.01A



Square
EFF2104.01B



D2.3 IMPLACIL® 5.0

INTERNAL CONNECTION




Ucla Plastic Surgery with CrCo



Hexagonal
EFF2104.01A



Square
EFF2104.01B

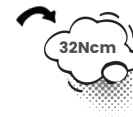

Link to CAD/CAM TiBase



Hexagonal
EFF2104.01A



Square
EFF2104.01B



Interface



Hexagonal
EFF2104.03A





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal!
Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2102.14A (O2.2)
(ø1.6mm thread)
EFF2103.15A (B2.1)
(ø1.8mm thread)



Square EFF2102.14B (O2.2)
(ø1.6mm thread)
EFF2103.15B (B2.1)
(ø1.8mm thread)



O2.2
FRIALIT®
3.8



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2102.14A (O2.2)
(ø1.6mm thread)
EFF2103.15A (B2.1)
(ø1.8mm thread)



Square EFF2102.14B (O2.2)
(ø1.6mm thread)
EFF2103.15B (B2.1)
(ø1.8mm thread)



B2.1
SIN®
3.8



Link to CAD/CAM TiBase



Hexagonal
EFF2102.14A (O2.2)
(ø1.6mm thread)
EFF2103.15A (B2.1)
(ø1.8mm thread)



Square EFF2102.14B (O2.2)
(ø1.6mm thread)
EFF2103.15B (B2.1)
(ø1.8mm thread)



INTERNAL CONNECTION



Interface



Hexagonal EFF2102.15A (O2.2)
(ø1.6mm thread)
EFF2103.11A (B2.1)
(ø1.8mm thread)





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2102.14A (O2.3)
(ø1.6mm thread)
EFF2103.15A (B2.2)
(ø1.8mm thread)



Square EFF2102.14B (O2.3)
(ø1.6mm thread)
EFF2103.15B (B2.2)
(ø1.8mm thread)



02.3
FRIALIT®
45



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2102.14A (O2.3)
(ø1.6mm thread)
EFF2103.15A (B2.2)
(ø1.8mm thread)



Square EFF2102.14B (O2.3)
(ø1.6mm thread)
EFF2103.15B (B2.2)
(ø1.8mm thread)



B2.2
SIN®
45



Link to CAD/CAM TiBase



Hexagonal
EFF2102.14A (O2.3)
(ø1.6mm thread)
EFF2103.15A (B2.2)
(ø1.8mm thread)



Square EFF2102.14B (O2.3)
(ø1.6mm thread)
EFF2103.15B (B2.2)
(ø1.8mm thread)



INTERNAL CONNECTION



Interface



Hexagonal EFF2102.15A (O2.3)
(ø1.6mm thread)
EFF2103.11A (B2.2)
(ø1.8mm thread)





Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw




Temporary Abutment



Hexagonal
EFF2102.14A (O2.4)
(ø1.6mm thread)
EFF2103.15A (B2.3)
(ø1.8mm thread)



Square EFF2102.14B
(O2.4)
(ø1.6mm thread)
EFF2103.15B (B2.3)
(ø1.8mm thread)

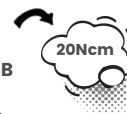

O2.4
FRIALIT®
5.5



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2102.14A (O2.4)
(ø1.6mm thread)
EFF2103.15A (B2.3)
(ø1.8mm thread)



Square EFF2102.14B
(O2.4)
(ø1.6mm thread)
EFF2103.15B (B2.3)
(ø1.8mm thread)



B2.3
SIN®
5.5



Link to CAD/CAM TiBase



Hexagonal
EFF2102.14A (O2.4)
(ø1.6mm thread)
EFF2103.15A (B2.3)
(ø1.8mm thread)



Square EFF2102.14B
(O2.4)
(ø1.6mm thread)
EFF2103.15B (B2.3)
(ø1.8mm thread)



INTERNAL CONNECTION





Interface



Hexagonal EFF2102.15A
(O2.4)
(ø1.6mm thread)
EFF2103.11A (B2.3)
(ø1.8mm thread)



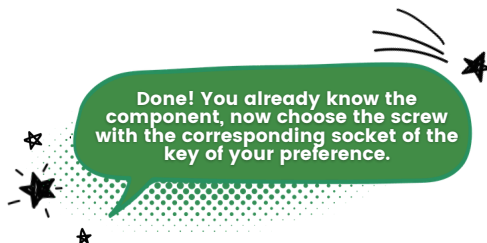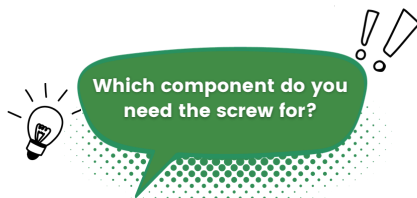

Hexagonal
EFF2104.01A



Square
EFF2104.01B



G2.1
DENTOFLEX®
325



Dynamic Link



Hexagonal
EFF2104.01A




Square
EFF2104.01B





INTERNAL CONNECTION



Hexagonal
EFF2104.03A

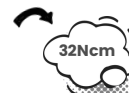


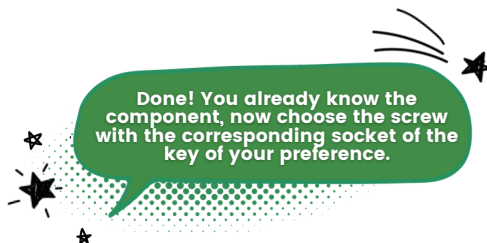
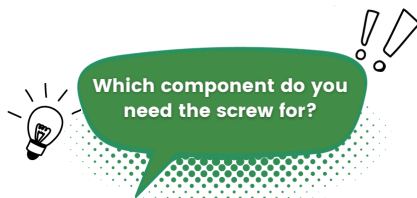


G2.2
DENTOFLEX®
3.5



INTERNAL CONNECTION

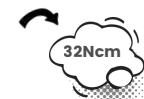

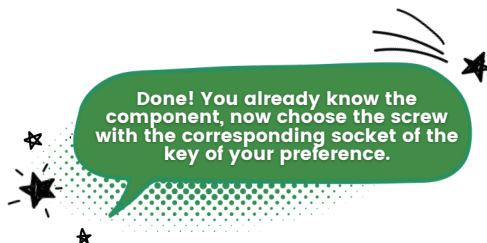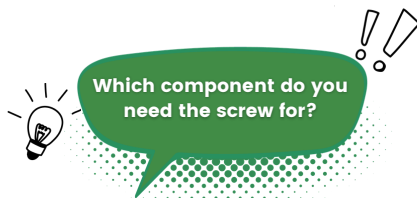

Hexagonal
EFF2104.01A



Square
EFF2104.01B



G2.3
DENTOFLEX®
3.75



Dynamic Link



Hexagonal
EFF2104.01A




Square
EFF2104.01B






INTERNAL CONNECTION



Hexagonal
EFF2104.03A

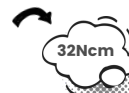

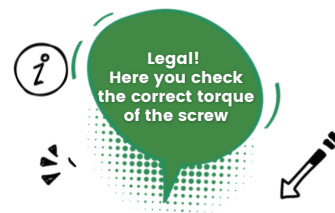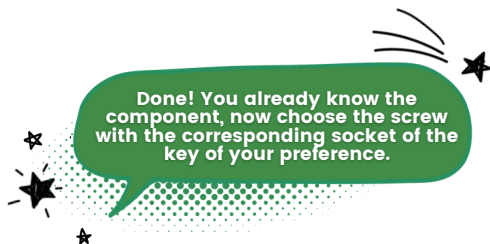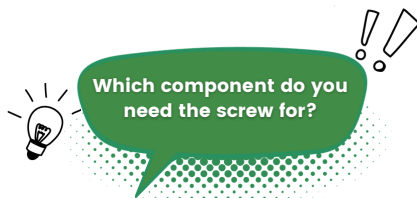

P2.1

CAMLOG®
3.3

INTERNAL
CONNECTION



Temporary
Abutment



Hexagonal
EFF2102.04A





Link to
CAD/CAM
TiBase



Hexagonal
EFF2102.04A



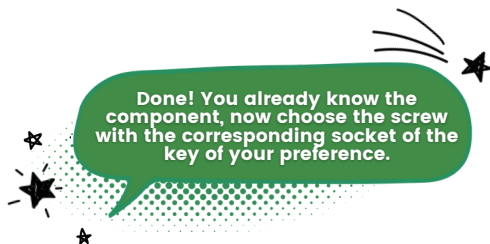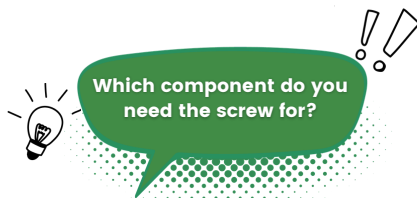

P2.2

CAMLOG®
3.8

INTERNAL CONNECTION



Temporary
Abutment



Hexagonal
EFF2104.06A






Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.06A



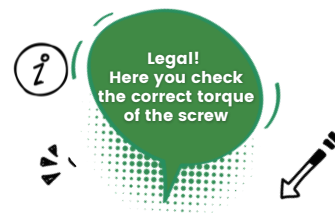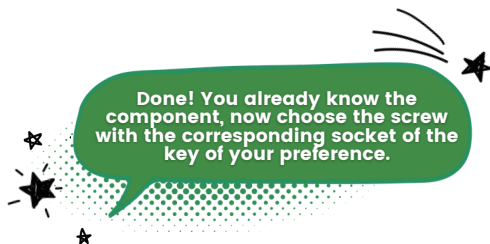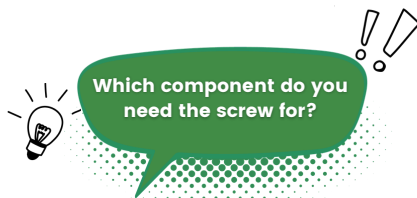

P2.3

CAMLOG®
4.3

INTERNAL CONNECTION



Temporary
Abutment



Hexagonal
EFF2104.06A



32Ncm






Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.06A



32Ncm

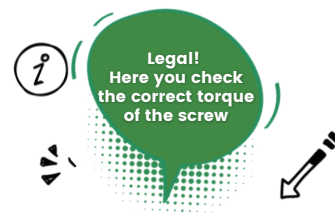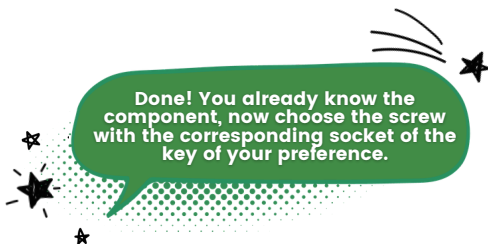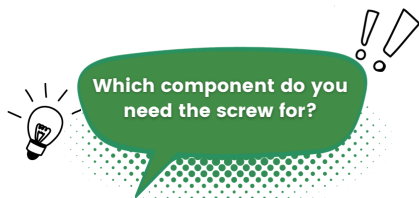

Sistema Multi Plataforma

P2.4

CAMLOG®
5.0

INTERNAL CONNECTION



Temporary
Abutment



Hexagonal
EFF2104.06A






Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.06A



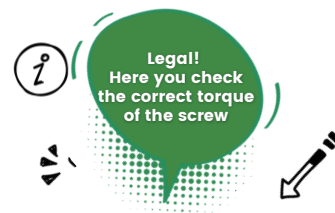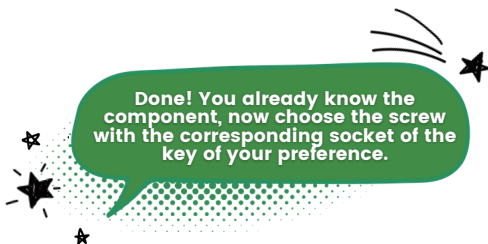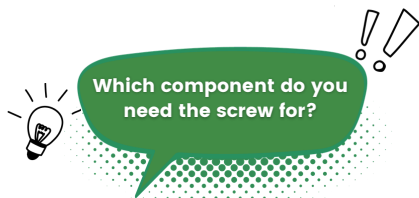

Temporary
Abutment



Hexagonal
EFF2104.06A




P2.5

CAMLOG®
6.0

INTERNAL CONNECTION





Link to
CAD/CAM
TiBase

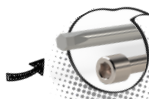

Hexagonal
EFF2104.06A



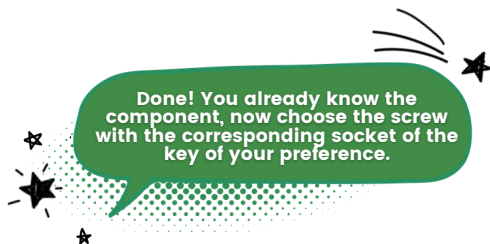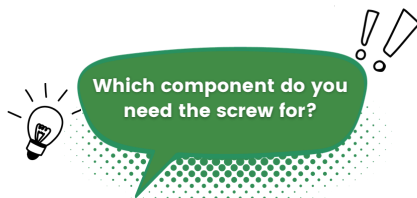

Sistema Multi Piattaforma

X2.1

BAUMER®
4.0

INTERNAL CONNECTION



Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.01A



Square
EFF2104.01B



Dynamic Link





Hexagonal
EFF2104.01A



Square
EFF2104.01B



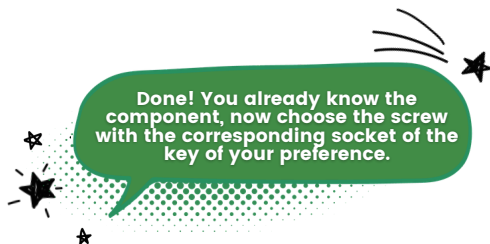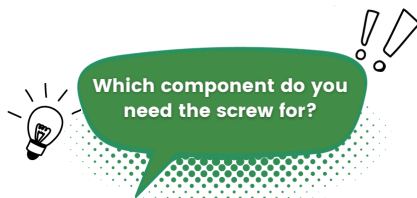

C2.1
C2.2
C2.3

CONNECTION®
4.1/4.3/5.0/6.0

INTERNAL
CONNECTION



Temporary
Abutment



Hexagonal
EFF2103.18A



Square EFF2103.18B



Link to
CAD/CAM
TiBase





Hexagonal
EFF2103.18A



Square EFF2103.18B



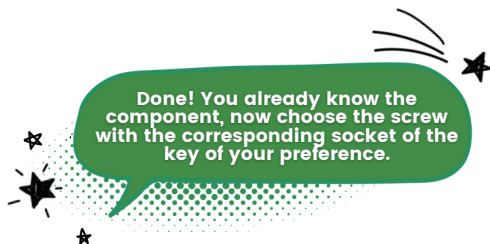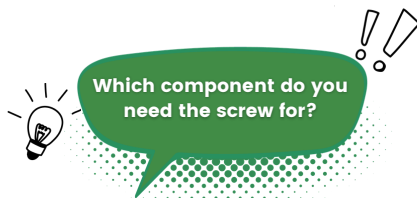

Link to
CAD/CAM
TiBase



Hexagonal
EFF2103.25A



E2.1

TITANIUMFIX®
3.75/3.0

INTERNAL
CONNECTION





Dynamic Link



Hexagonal
EFF2103.25A




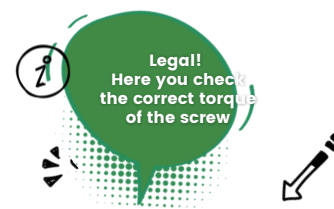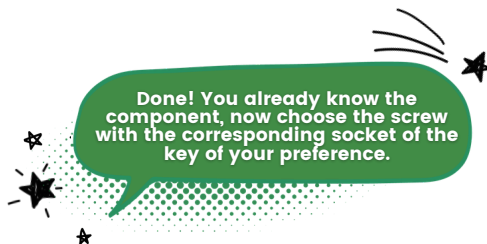

L2.1

BIOMET 3I™
3.4

INTERNAL CONNECTION

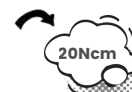



Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal! Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2102.06A



L2.2

BIOMET 3I™
4.1

INTERNAL CONNECTION



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2102.06A



Link to
CAD/CAM
TiBase



Hexagonal
EFF2102.06A





Interface



Hexagonal
EFF2102.15A



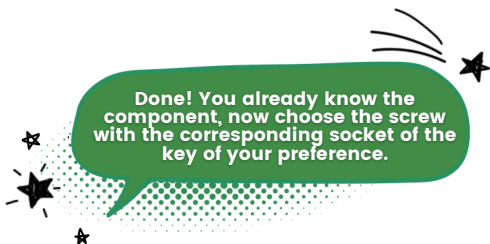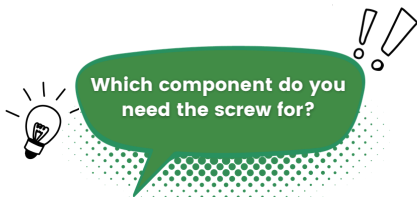


L2.3

BIOMET 3I™
5.0

INTERNAL CONNECTION



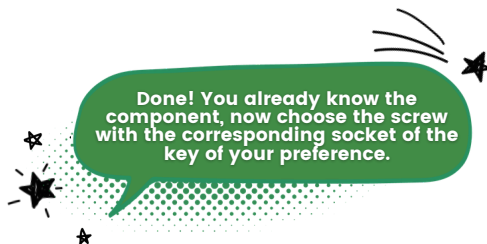

A2.1


NEODENT®
4.3|5.0

INTERNAL CONNECTION





Prosthetic Kit

All in One / All in One Full / Easy

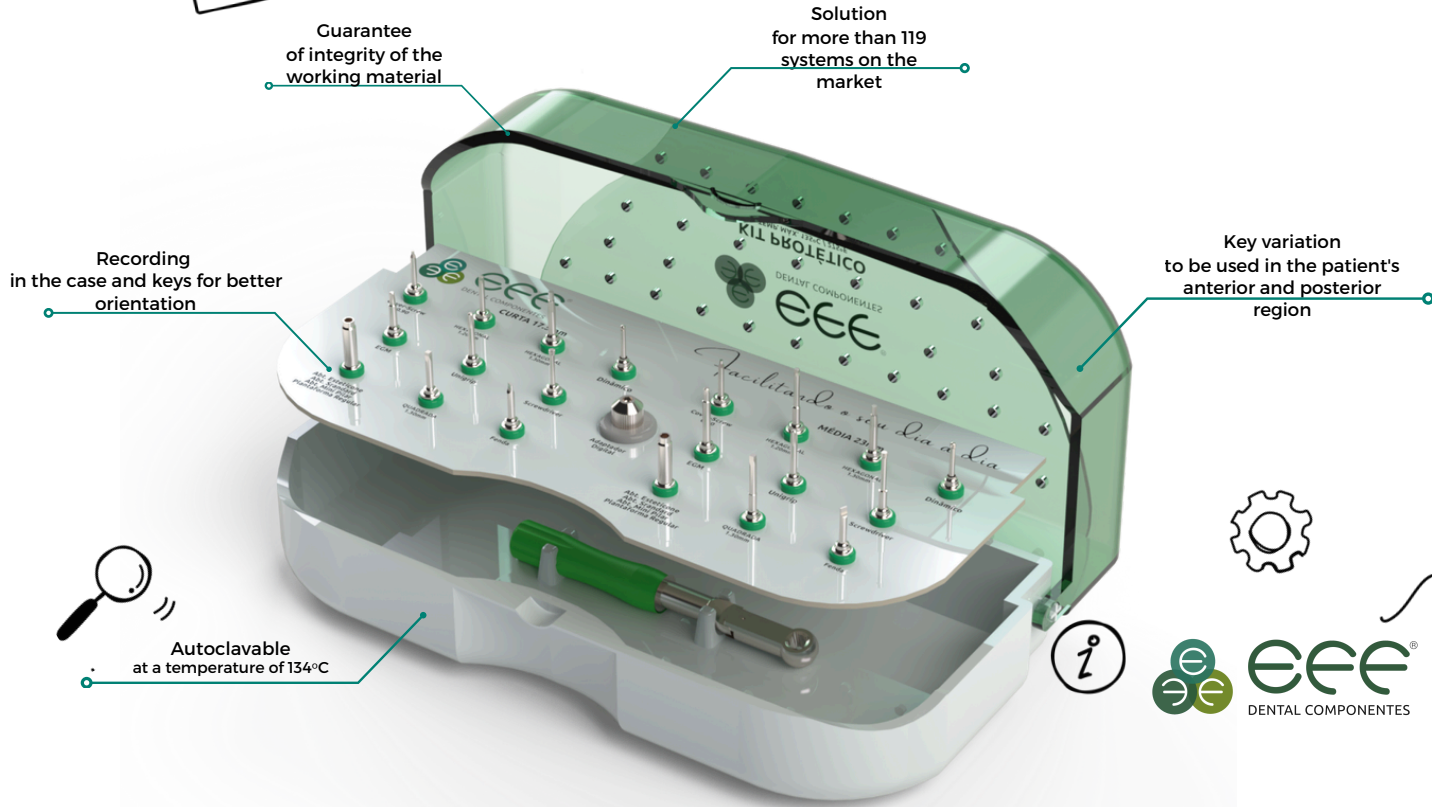

eee
DENTAL COMPONENTES



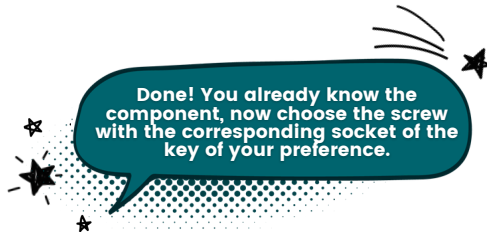
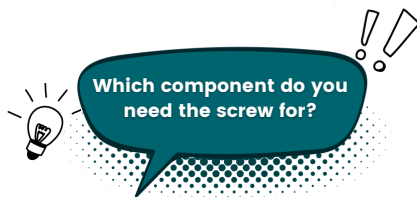


K1.1


BRANEMARK
SYSTEM® 3.5

EXTERNAL
CONNECTION



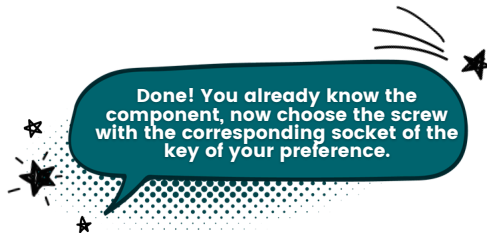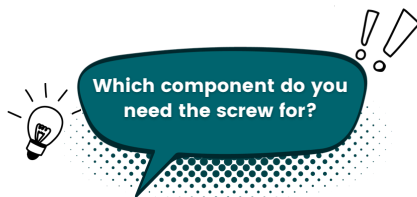

Temporary Abutment

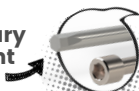

Hexagonal
EFF2104.01A



Grip
EFF2104.01C



Square
EFF2104.01B



Dynamic
EFF2104.01F



Systema Multi-Platform



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2104.01A



Grip
EFF2104.01C



Square
EFF2104.01B



Dynamic
EFF2104.01F



K1.2

BRANEMARK SYSTEM® 4.1

EXTERNAL CONNECTION



Prov. EVOLUTION



Hexagonal
EFF2104.01A



Grip
EFF2104.01C



Square
EFF2104.01B



Dynamic
EFF2104.01F



Ucla Plastic



Hexagonal
EFF2104.01A




Grip
EFF2104.01C



Square
EFF2104.01B



Dynamic
EFF2104.01F

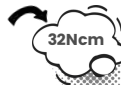

Link to CAD/CAM TiBase



Hexagonal
EFF2104.01A



Grip
EFF2104.01C



Square
EFF2104.01B



Dynamic
EFF2104.01F



Dynamic Link



Hexagonal
EFF2104.01A



Grip
EFF2104.01C





Square
EFF2104.01B



Dynamic
EFF2104.01F



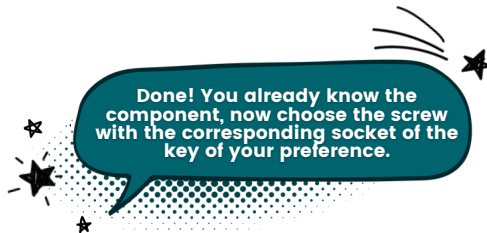
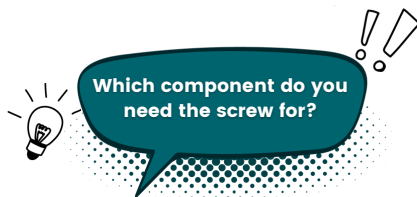


K1.3

BRANEMARK SYSTEM® 3.5

EXTERNAL CONNECTION



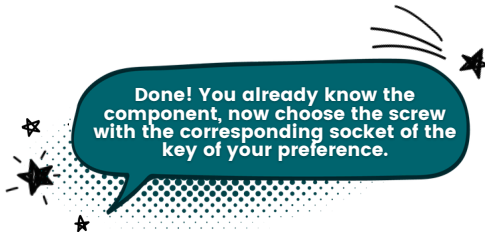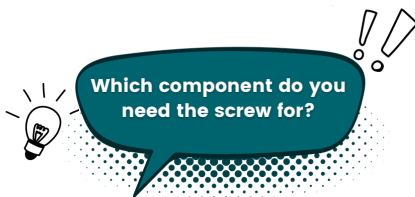

Temporary Abutment



Hexagonal EFF2103.02A



Ucla Plastic Surgery with CrCo



Hexagonal EFF2103.02A



A1.1

NEODENT® 3.3



Prov. EVOLUTION



Hexagonal EFF2103.02A



C1.1

CONNECTION® 3.3



Ucla Plastic



Hexagonal EFF2103.02A



EXTERNAL CONNECTION



Link to CAD/CAM TiBase



Hexagonal EFF2103.02A





Dynamic Link



Hexagonal EFF2103.02A



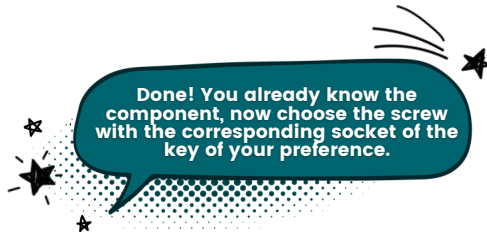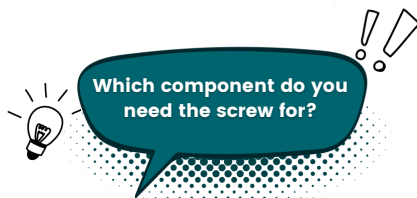

**Temporary
Abutment**



**Hexagonal
EFF2104.01A**



32Ncm



A1.3

**NEODENT®
4.3**

EXTERNAL CONNECTION



**Ucla Plastic
Surgery with CrCo**



**Hexagonal
EFF2104.01A**



32Ncm



Ucla Plastic



**Hexagonal
EFF2104.01A**



32Ncm



**Link to
CAD/CAM
TiBase**



**Hexagonal
EFF2104.01A**



32Ncm



Which component do you need the screw for?



Done! You already know the component, now choose the screw with the corresponding socket of the key of your preference.



Legal!
Here you check the correct torque of the screw



Temporary Abutment



Hexagonal
EFF2105.02A



Square
EFF2105.02B



32Ncm

A1.4
NEODENT®
5.0



Ucla Plastic Surgery with CrCo



Hexagonal
EFF2105.02A



Square
EFF2105.02B



32Ncm

C1.3
CONNECTION®
5.0
EXTERNAL CONNECTION



Ucla Plastic



Hexagonal
EFF2105.02A



Square
EFF2105.02B



32Ncm



Link to CAD/CAM TiBase





Hexagonal
EFF2105.02A



Square
EFF2105.02B



32Ncm

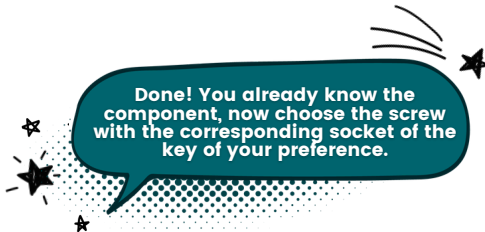
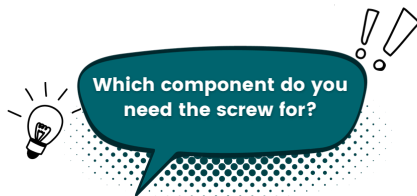


L1.1
BIOMET 3I™
3.4



B1.1
SIN®
3.4
EXTERNAL CONNECTION



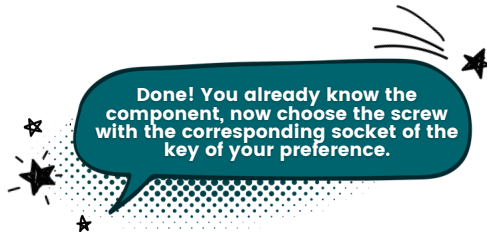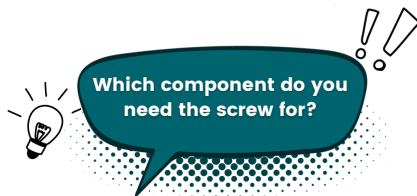

Temporary Abutment



Hexagonal EFF2104.01A



L1.3

BIOMET 3I™ 5.0



Ucla Plastic Surgery with CrCo



Hexagonal EFF2104.01A



32Ncm

B1.3

SIN® 5.0

EXTERNAL CONNECTION



Ucla Plastic



Hexagonal EFF2104.01A



32Ncm





Link to CAD/CAM TiBase



Hexagonal EFF2104.01A



32Ncm

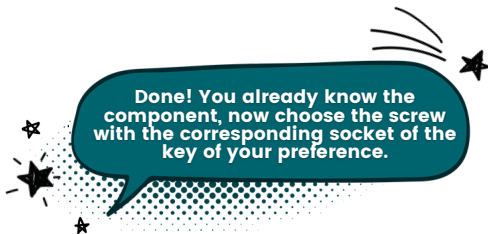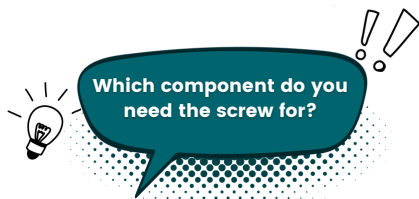

D1.1

IMPLACIL®
3.3

EXTERNAL
CONNECTION



Temporary
Abutment



Hexagonal
EFF2103.02A



Square
EFF2103.02BB



Link to
CAD/CAM
TiBase





Hexagonal
EFF2103.02A



Square
EFF2103.02BB



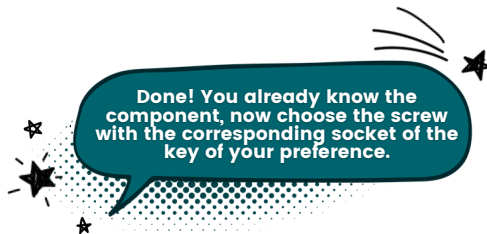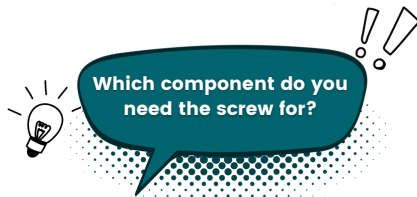

Sistema Multi Piattaforma

E1.1

TITANIUMFIX® 3.3

EXTERNAL CONNECTION



**Temporary
Abutment**



**Hexagonal
EFF2106.01A**





**Link to
CAD/CAM
TiBase**



**Hexagonal
EFF2106.01A**




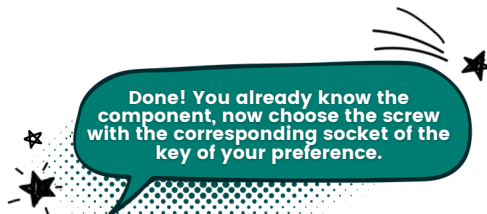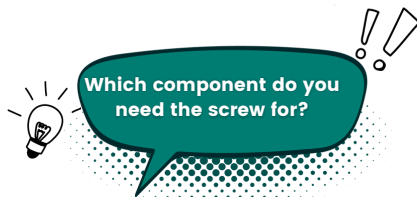

Plastic Coping



Hexagonal EFF2101.01A
(h 4.1mm)
EFF2101.01A.SAE (h 4.6mm)

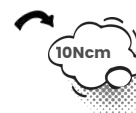

CrCo Base
Plastic Coping




Hexagonal EFF2101.01A
(h 4.1mm)
EFF2101.01A.SAE (h 4.6mm)



Link to
CAD/CAM
TiBase





Hexagonal EFF2101.01A
(h 4.1mm)
EFF2101.01A.SAE (h 4.6mm)

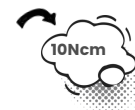

M.P.
MULTI-UNIT -
REGULAR
PLATFORM



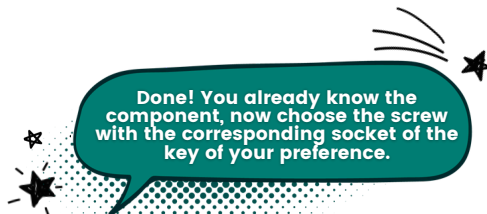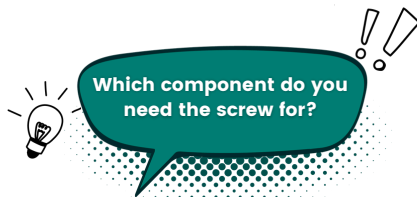

Plastic Copping



Hexagonal
EFF2103.01A



Grip EFF2103.01C



M.P.L.1

MULTI-UNIT
PLATFORM
LONG



CrCo Base
Plastic Copping



Hexagonal
EFF2103.01A



Grip EFF2103.01C



Link to
CAD/CAM
TiBase





Hexagonal
EFF2103.01A



Grip EFF2103.01C




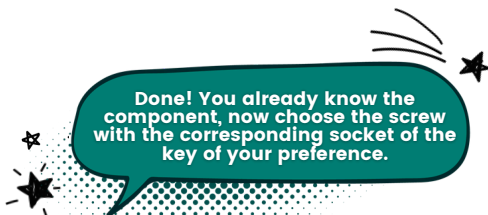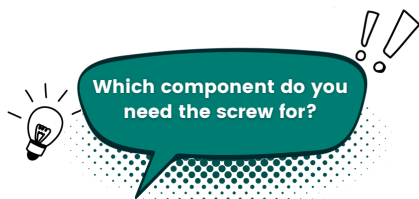

Plastic Coping




Hexagonal
EFF2104.22A

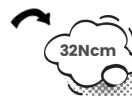

CrCo Base
Plastic Coping




Hexagonal
EFF2104.22A

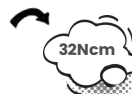

Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.22A

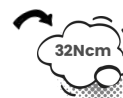

M.P.L.2

MULTI-UNIT
PLATFORM
LONG - NEO






CrCo Base
Plastic Coping



Hexagonal
EFF2104.38A



Dynamic
EFF2104.38F

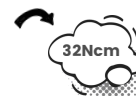

Link to
CAD/CAM



Hexagonal
EFF2104.38A



Dynamic
EFF2104.38F




Link to
CAD/CAM
TiBase



Hexagonal
EFF2104.38A





Dynamic
EFF2104.38F

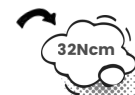

EONE
AESTHETICS




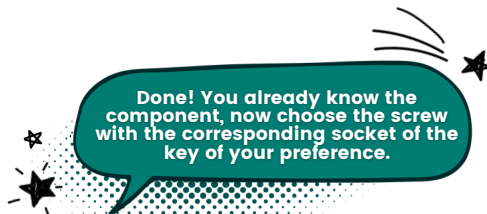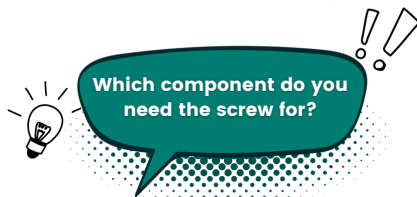

Plastic Coping




Hexagonal
EFF2101.01A.SAE
(pair h4.6mm)

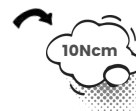

CrCo Base
Plastic Coping



Hexagonal
EFF2101.01A.SAE
(pair h4.6mm)

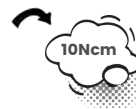

Link to
CAD/CAM
TiBase



Hexagonal
EFF2101.01A.SAE
(pair h4.6mm)



EAST
ESTHETICONE





**PLEASE CONTACT US FOR MORE
INFORMATION ABOUT SCREWS,
AT**

 (11) 3059-2200

 www.efdental.com.br

 @efdental

MKT8.1 EFF DENTAL_EFF DENTAL SCREW SELECTION MANUAL_Rev01_DATE: 15/09/2023

All rights reserved. EFF Dental, the EFF Dental logo and all other trademarks used herein are, unless otherwise stated or if not evident from the context of the case, trademarks of EFF Dental. Product images are not necessarily to scale. Some products may not be regulatory approved/cleared in all markets. Please contact your local EFF Dental sales office for information on product availability.